

SALTER LABS®

Salter-Style® **Demand Cannulas**

Salter-Style® cannulas specifically designed for use with Dual Port Oxygen Conservation Delivery systems. Adult, Small Adult/Pediatric sizes with 4', 5', 7' supply tubes.



Features

- ☐ Dual tubing set
- ☐ Soft, curved nasal prongs
- ☐ Anatomically curved facepiece to fit upper lip
- ☐ Over-the-ear style
- ☐ *TLCannula™* model has soft foam cushions on headset tubes
- ☐ Available in Small Adult/Pediatric and Adult sizes with 4', 5', 7' tube lengths

Benefits

- ☐ Provides separate independent tubes for oxygen delivery and sensing inspiratory/expiratory effort
- ☐ Improved comfort and anatomical fit
- ☐ Provides a "natural" fit, lowers patient anxiety
- ☐ A secure fit, easy to put on or take off, even for a reclining patient
- ☐ The foam cushions provide extra comfort for the over-the-ear area and aid in holding the cannula securely in place
- ☐ Select the appropriate model for your stationary or ambulatory patients

**Variety of uses. . . Hospitals, Skilled Nursing Facility, Clinics, Homecare,
Assisted Living, Extended Care**




slc 08

*Salter-Style® Demand Cannulas**

Special cannulas designed to give both adult and pediatric users of dual port oxygen conservation delivery systems all the comfort and conveniences of a *Salter-Style®* cannula. Unique dual tubing delivers oxygen on demand through one tube while inspiratory / expiratory effort is sensed through the other tube.

Sense & Delivery to Separate Nostrils


Permanent Barrier

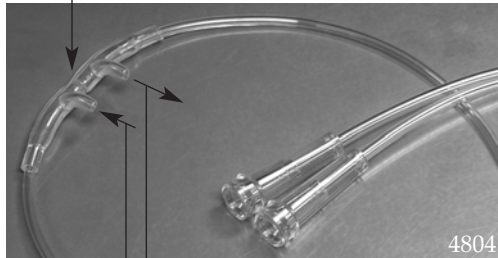

Patient Sensing *Oxygen Delivery*

Sense & Delivery to Both Nostrils

Patient Sensing and Oxygen Delivery

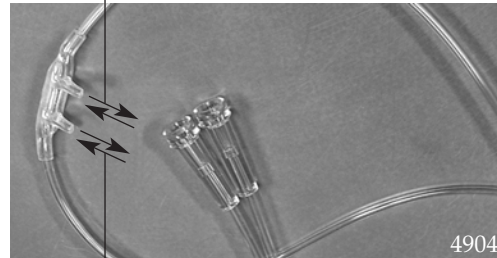

Patient Sensing and Oxygen Delivery

The *Salter-Style® DEMAND TLCannula™*. Setting the standard in patient comfort and convenience for long term oxygen therapy (LTOT) patients.

The Demand *TLCannula™* is assembled with gray colored, closed cell, soft foam cushions placed over the cannula head tubes. The cushions provide a new dimension in comfort for sensitive over-the-ear tissue.

The closed cell foam cushions are also available as separate items for self installation on nasal cannulas. (See SLC-71 Salter Accessories E-Z Wrap™ P/N 1016)

Ordering Information

*Salter-Style® Demand Cannula**

*Designed for dual port oxygen conservation delivery systems.
For flow rates up to 6 LPM*



Gray colored, closed cell, soft foam cushions are pre assembled over the headset tubing.

TLCannula™ P/N 4804TLC

	Units (per case)	Part Number
<i>Sense and Delivery to separate nostrils</i>		
Adult , Demand cannula with 4' supply tubes	25	4804
Adult , Demand cannula with 5' supply tubes	25	4805
Adult , Demand cannula with 7' supply tubes	25	4807
Small Adult/Pediatric Demand cannula with 4' supply tubes	25	4824
Small Adult/Pediatric Demand cannula with 7' supply tubes	25	4827
Adult <i>Salter-Style®</i> Demand <i>TLCannula™</i> with foam cushions and 4' supply tubes	10 or 25	4804TLC
Adult <i>Salter-Style®</i> Demand <i>TLCannula™</i> with foam cushions and 7' supply tubes	10 or 25	4807TLC
<i>Sense and Delivery to both nostrils</i>		
Adult , Dual Delivery / Dual Sense Demand cannula, 4' tubes	25	4904
Adult , Dual Delivery / Dual Sense Demand cannula, 5' tubes	25	4905
Adult , Dual Delivery / Dual Sense Demand cannula, 7' tubes	25	4907