


SALTER LABS®

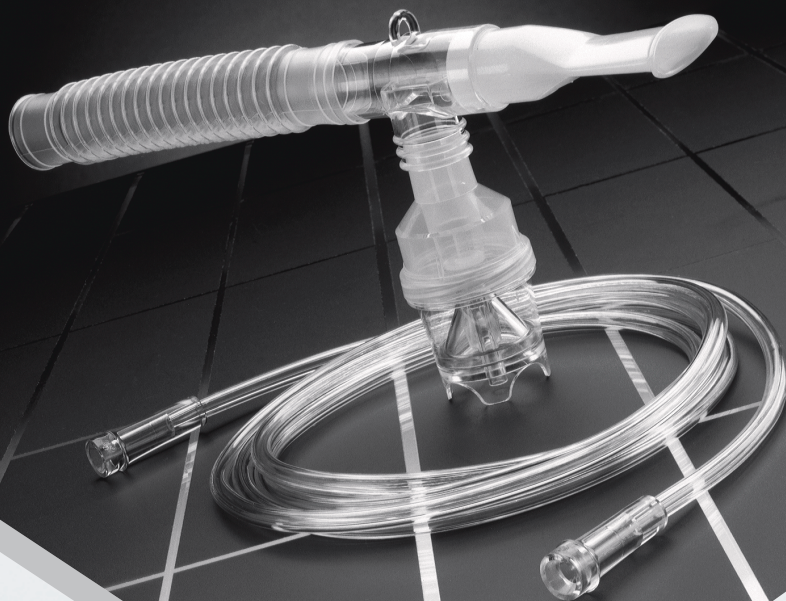

Salter 8900 Series

Small Volume Jet Nebulizer

Disposable—single patient use

For reliability, efficiency and convenient aerosol treatment.



Features

- ☐ Nebulizes 3cc in less than 10 minutes
- ☐ Nebulizes in horizontal or vertical position
- ☐ Larger surface area provided by unique convex cone design
- ☐ Removable green convex cone
- ☐ Scalloped, bottom cuff design
- ☐ One-eighth turn screw on cap
- ☐ Anti-drool "T"

Benefits

- ☐ Shortened treatment times
- ☐ Allows patient to be in a comfortable position for treatment
- ☐ Maximum capillary action reduces treatment times
- ☐ Easy to see, and allows for thorough jet cleaning and reassembly
- ☐ Allows the unit to stand alone; simplifies tube attachment
- ☐ Easy to assemble and disassemble, provides a leak proof seal
- ☐ Prevents medication contamination from patient saliva

Variety of uses. . . Hospital, Skilled Nursing Facility, Homecare, Clinic




slc-21

Salter 8900 Series Jet Nebulizer*

Disposable—single patient use

Built-in quality and reliability,
minimizing patient operational
assistance.

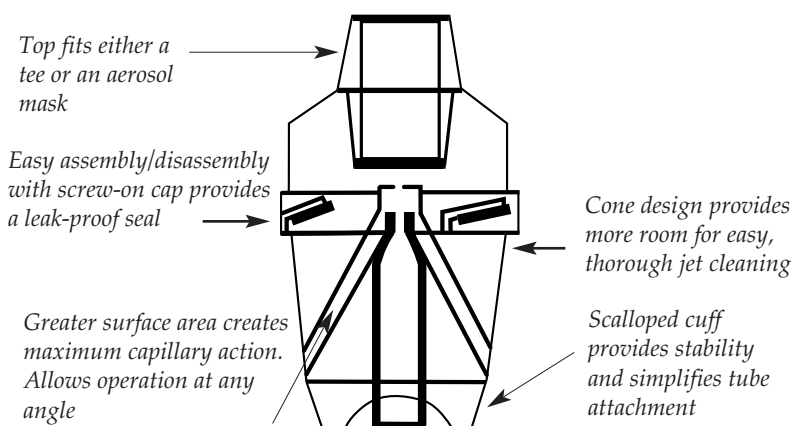

Flow Rate (L/min.)	Particle Size (microns)										MMAD
	0.14	0.23	0.40	0.62	1.30	2.30	3.20	4.60	6.50	9.20	
4	3%	5%	10%	21%	25%	26%	6%	2%	1%	2%	1.30
6	4%	9%	15%	26%	25%	16%	2%	1%	2%	1%	1.10
7	4%	11%	16%	27%	22%	14%	2%	1%	1%	1%	1.10
8	4%	12%	17%	30%	20%	12%	2%	0%	1%	1%	1.00
10	4%	14%	19%	31%	17%	10%	1%	1%	1%	1%	0.98

Test Specifications: were performed at 50psi source pressure with 0.9%NaCl. Particle size distribution measured by QCM Cascade impactor.

Ordering Information

8900 Series Jet Nebulizer Kits*

Disposable—single patient use

	Units (per case)	Part No.		Units (per case)	Part No.
Full Kit: Nebulizer, anti-drool "T", mouthpiece, 6" reservoir tube, 7' supply tube	50	8900	Aerosol "T" 22mm O.D. x 22mm O.D. x 15mm I.D.	50	8917
Same as above with 7' supply tube with female thread-grip connector	50	8900TG	Nebulizer, adult over-the-ear style aerosol mask	50	8903
In-Line Kit: Nebulizer, 22mm O.D. x 22mm O.D. x 18mm I.D. "T", 6" reservoir tube, 7' supply tube	50	8915	Nebulizer, adult elastic strap style aerosol mask	50	8923
Same as above with 7' supply tube with female thread-grip connector	50	8915TG	Nebulizer, adult over-the-ear style aerosol mask, 7' supply tube	50	8904
In-Line "T" 22mm O.D. x 22mm O.D. x 18mm I.D.	50	8916	Same as above with female thread grip connector	50	8904TG
Nebulizer, anti-drool "T", mouthpiece, 7' supply tube	50	8911	Nebulizer, adult elastic strap style aerosol mask, 7' supply tube	50	8924
Nebulizer, anti-drool "T", mouthpiece, 6" reservoir tube	50	8912	Nebulizer, pediatric elastic strap style aerosol mask	50	8905
Nebulizer, anti-drool "T", mouthpiece	50	8913	Nebulizer, pediatric elastic strap style aerosol mask, 7' supply tube	50	8906
Nebulizer, anti-drool "T", 6" reservoir tube	50	8914	Nebulizer, pediatric aerosol mask VADS*, with elastic head strap, and 7' supply tube	50	8907
Nebulizer, 7' supply tube	50	8902			
Nebulizer only	50	8901			
Anti-Drool "T" 22mm O.D. x 22mm I.D. x 18mm I.D.	50	8910			
Mouthpiece	50	1003			

*VADS (Valved Aerosol Delivery System SLC-43)

Printed in U.S.A. *U.S. Patent No. RE33,717 © Copyright, 2009 Salter Labs®
Revised January, 2011 SLC-21

SALTER LABS®

100 W. Sycamore Road - Arvin, CA 93203
Phone: 661-854-3166
1-800-421-0024
1-800-235-4203
Fax: 661-854-3850
Toll Free Fax: 1-800-628-4690
Fax: 661-854-3850
E-mail: salterlabs@us.salterlabs.com
www.salterlabs.com

CE
0482

ISO 13485