

SALTER LABS®



Humidifier

*P/N 7100 and P/N 7600 350 cc
Bubble Humidifier - Dry*

360° Rectangular Diffuser ports

*provide quiet operation with low system back pressure.
Dome lid design combines high performance with rugged
durability. All Salter bubble humidifiers provide, dependable,
cost-effective performance. Ideal for use with all oxygen sources:
concentrators, liquid, tanks, wall outlets.*

Available with 3 PSI or 6 PSI Pop-off Alarm.



Features

- ☐ 360° molded, low resistance, high output diffuser head design
- ☐ Connecting tower design
- ☐ Large, quad wing nut connector
- ☐ Audible pop-off alarm
3PSI or 6PSI
- ☐ Jar permanently etched, large print minimum/maximum water level markings
- ☐ Easy seal, translucent, double helix thread jar

Benefits

- ☐ Quiet, efficient hydration with less rain out.
Resists clogging due to mineral build up from water.
- ☐ Greater resistance to damage from incidental impact
- ☐ Easy to grip, reduces chance of cross-threading
- ☐ Reduces the likelihood of concentrator damage.
Audible warning of downstream occlusion.
- ☐ Easy to read, durable will not rub, wear or wash off
- ☐ Simple, sure, half turn on/off with positive seal.
No gasket is required


Variety of uses... Homecare, Transport, Hospital, Extended Care Facilities, Hospice



slc-31A


Salter Humidifiers P/N 7100 3PSI pop off P/N 7600 6PSI pop off

Disposable Humidifiers from Salter Labs® are **single patient use products designed and engineered for long term durability and complete patient satisfaction.**

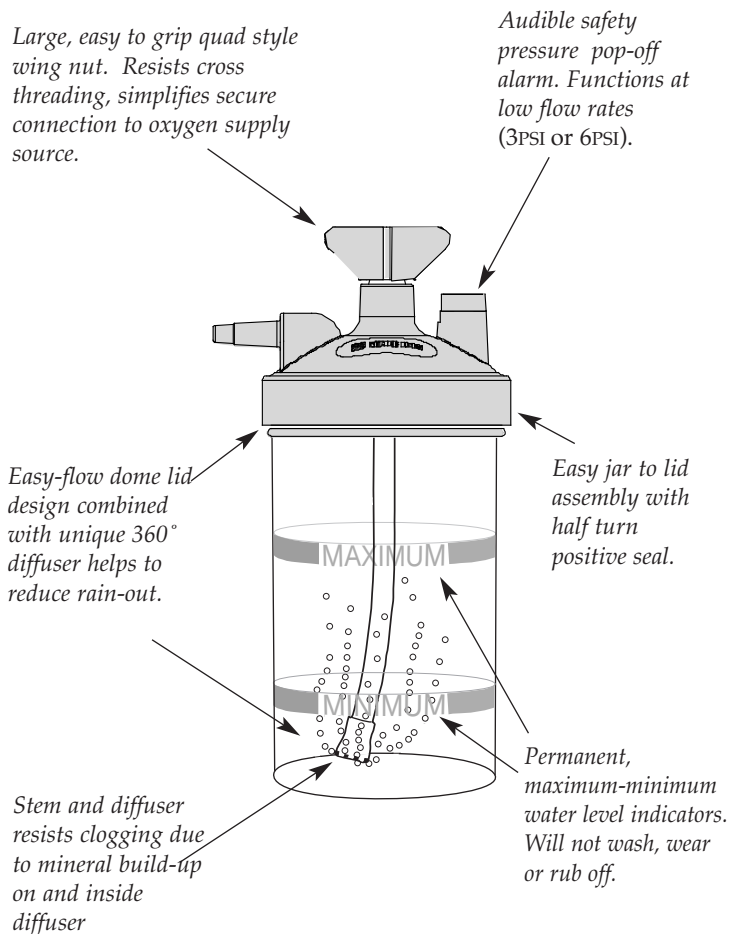

Start-up Kit

Available with P/N 7100 or P/N 7600 Humidifier

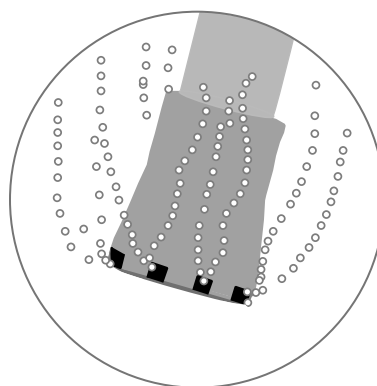

Unique 360° rectangular diffusion ports operate quietly with improved hydration and less system back pressure.

Ordering Information

Humidifier	Units (per case)		Part No.	
		3PSI	6PSI	
350 cc bubble humidifiers 3PSI (brown lid)	50	7100		
350 cc bubble humidifiers 6PSI (black lid)	50		7600	

For high flows (6-15 lpm) refer to our model #7900 Bubble Humidifier.

Start-up Kits

2 #1606 nasal cannulas, 1 #7100 or #7600 humidifier and 1 #2002 supply tube:

Kit with a 7' tube	20	7107	7607
Kit with a 25' supply tube	20	7125	7625
Kit with a 50' supply tube	20	7150	7650

Start-up Kits with smooth bore tubing available on request



100 W. Sycamore Road - Arvin, CA 93203
Phone: 661-854-3166
1-800-421-0024
1-800-235-4203
Fax: 661-854-3850
Toll Free Fax: 1-800-628-4690
Fax: 661-854-3850
E-mail: salterlabs@us.salterlabs.com
www.salterlabs.com

CE
0482